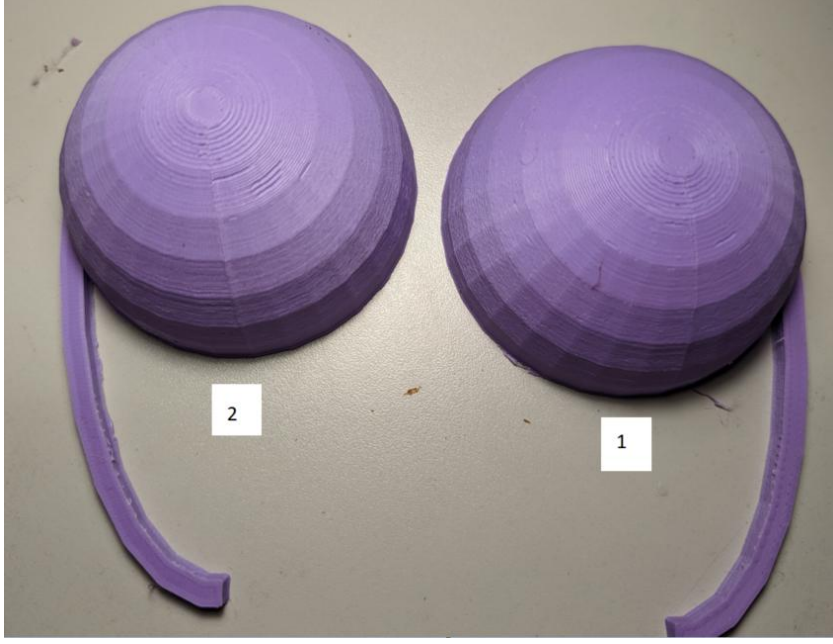


# Arthritis Typing Aid

## ASSEMBLY GUIDE

### Required Components

	<p><b>BOM</b></p> <ol style="list-style-type: none"><li>1. Aid-top</li><li>2. Aid-bottom</li></ol>
---	--

### Required Tools

- Superglue
- Sandpaper
- Soldering Iron (Optional)

### Required Personal Protective Equipment (PPE)

- Eye protection

# Arthritis Typing Aid

## ASSEMBLY GUIDE

### Assembly Instructions

#### Step 1

Take sandpaper and gently rough the flat bottom sides of the prints, then use super glue to glue them together as shown.

Alternatively, align the two halves of the device and run a soldering iron around the seam shown with the red dotted line to 'weld' the two halves together. This should be done in a well-ventilated room.

